


A INTUITIVE NEW PATENDED DESIGN


Only 2 easy steps are needed
for a secure assembly

ULTRA SAFETY PLUS TWIST

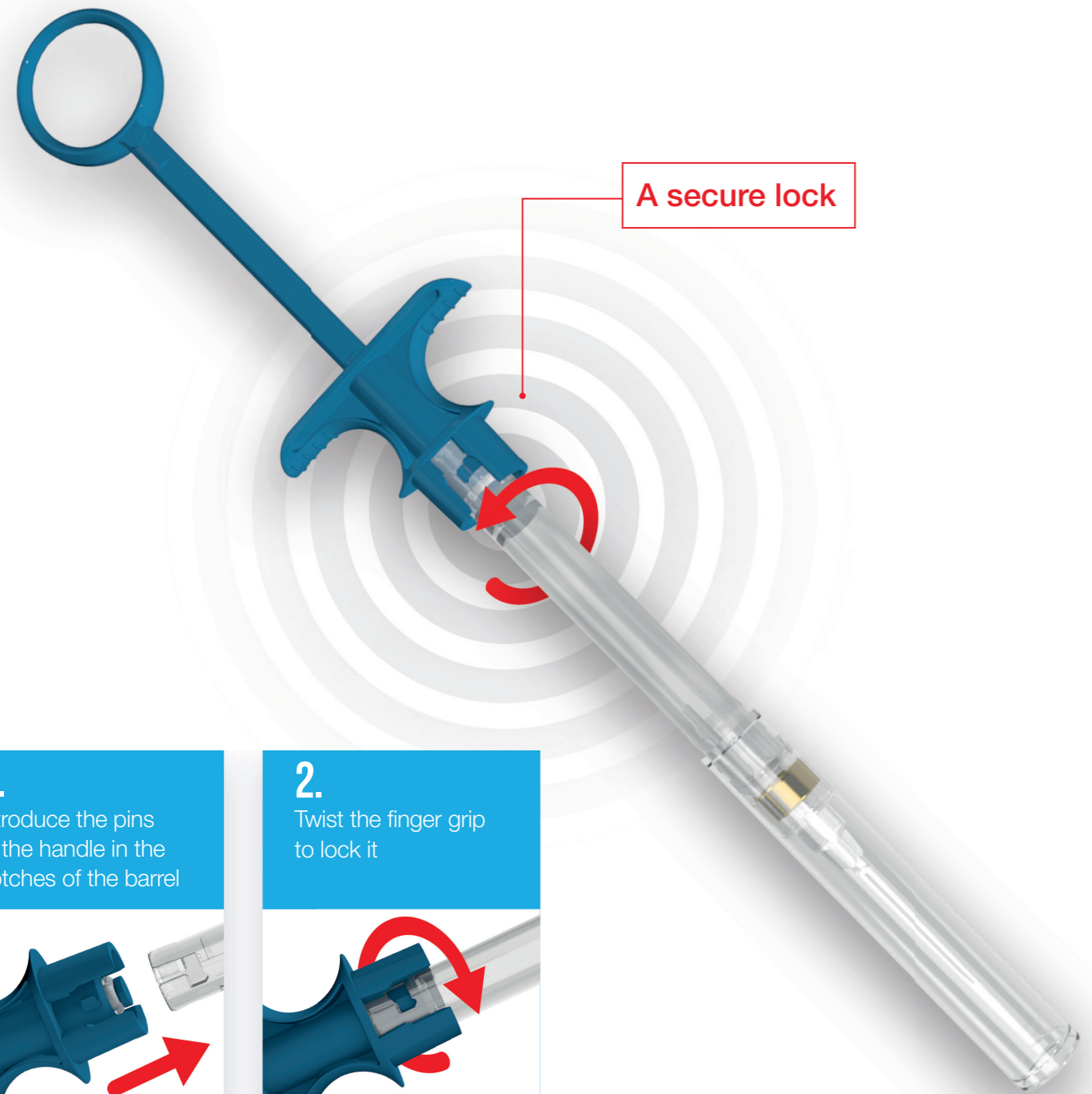


A secure lock

1.
Introduce the pins
of the handle in the
notches of the barrel

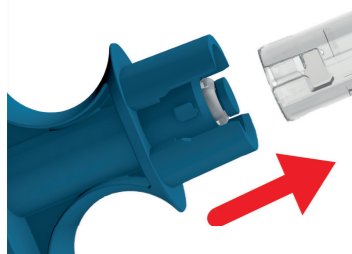

2.
Twist the finger grip
to lock it

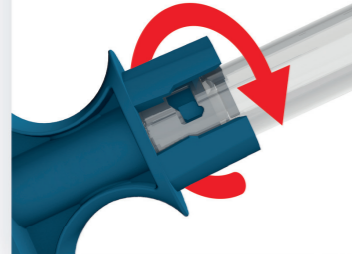

Find the relevant details in the instructions manual

ULTRA SAFETY PLUS TWIST



1 BOX OF 100 STERILE NEEDLES
+ 1 STERILISABLE HANDLE



1 BOX OF 50 STERILE SINGLE-USE HANDLES



1 BOX OF 50 STERILE NEEDLES



1 STERILISABLE HANDLE

LENGTH	GAUGE	COLOR CODE
X SHORT - 10 MM	30 G	
SHORT - 25 MM	30 G	
SHORT - 25 MM	27 G	
LONG - 35 MM	27 G	
LONG - 35 MM	25 G	

Please visit our website for more information:
www.septodont.com



ULTRA SAFETY PLUS TWIST

NEW!



TWIST
AND YOU'RE
COVERED



MANAGING
PAIN FOR
YOUR
PRACTICE



NEEDLE STICK INJURIES: A SIGNIFICANT RISK

Among the 35 million health workers, the World Health Organization estimated that 3 million suffered from percutaneous exposures to bloodborne pathogens each year*.



HEALTH

Among health care workers, about 40% of HBV and HCV infections and 2.5% of HIV infections are due to occupational sharps exposures*



COST

Costs for a healthcare organization when a needle stick injury occurs are about \$747 (range between \$199 and \$1,691)**



LEGISLATION


Many health agencies put in place legislations to prevent needle stick injuries requiring the implementation of safety procedures, the use of safety devices and the prohibition of the recapping of needles

ULTRA SAFETY PLUS

TWIST

SAFE & EASY

- > **Protects** you and your staff from needle stick injuries
- > **Complies** with latest regulations
- > **Intuitive** device
- > Available with either **sterile single use or sterilisable** handle

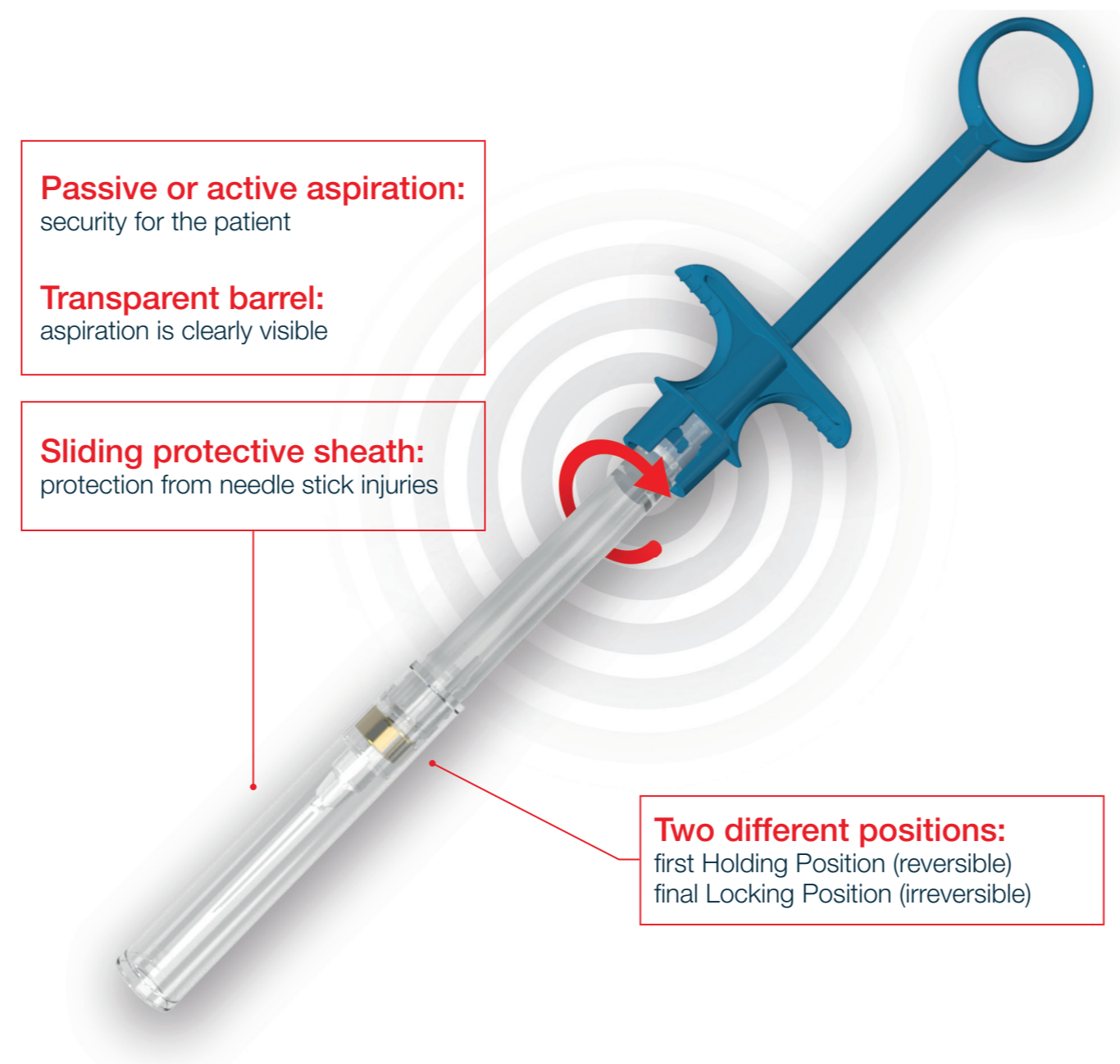

Passive or active aspiration:
security for the patient

Transparent barrel:
aspiration is clearly visible


Sliding protective sheath:
protection from needle stick injuries

Two different positions:
first Holding Position (reversible)
final Locking Position (irreversible)

A LONG-PROVEN EFFICACY

A study demonstrated that when introduced in a dental school, Ultra Safety Plus was a the key success factor for avoiding needle stick injuries.

FREQUENCY OF NEEDLE STICK INJURIES

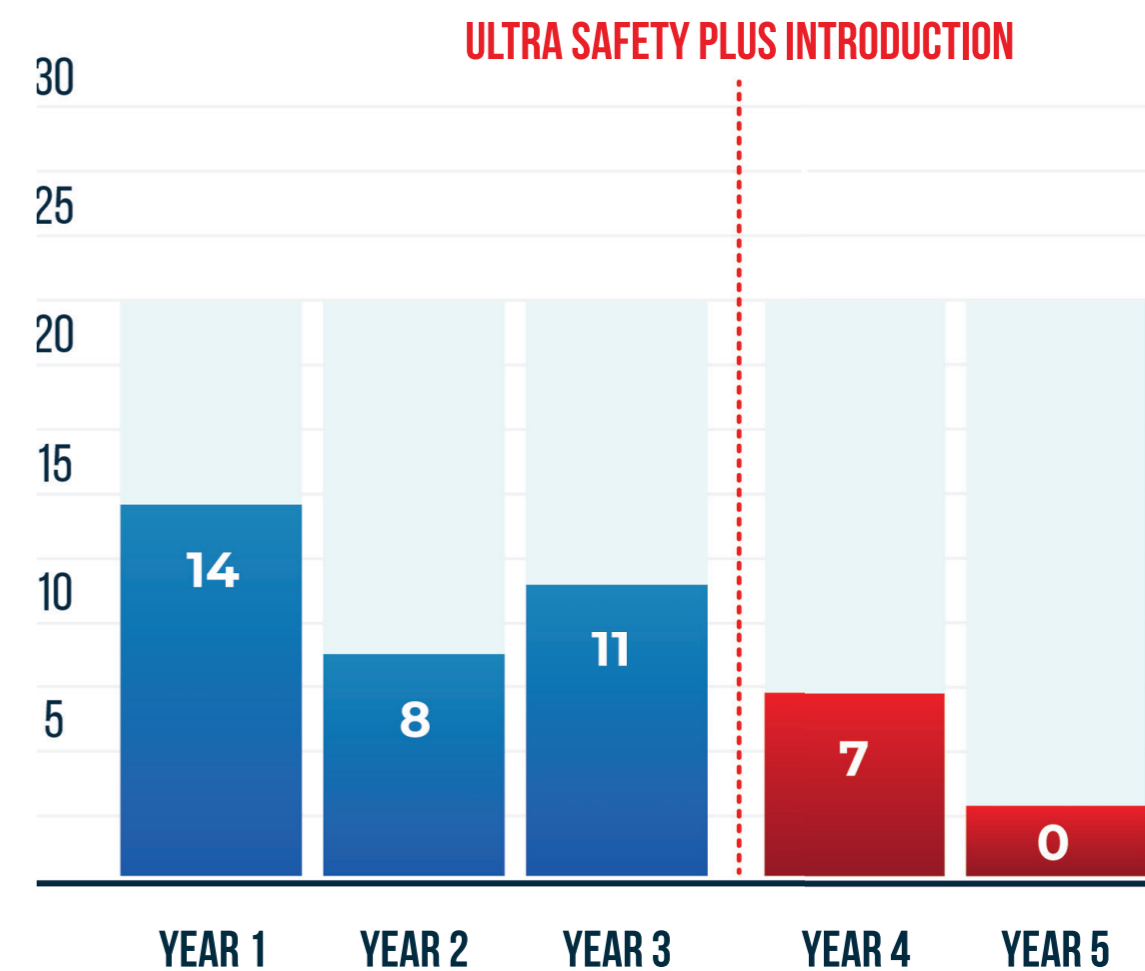

With Ultra Safety Plus, needle stick injuries decreased from an average of 11.8 to **0 injuries** per 1,000,000 hours worked⁽¹⁾

(1) J.M Zakrzewska et al. Introducing safety syringes into a UK dental school – a controlled study. Brit Dent J 2001 ; 190: 88-92.

* WHO, The world health report 2002 - Reducing Risks, Promoting Healthy Life, 2002
** Mannocci A, De Carli G, Di Bari V, et al. How Much do Needlestick Injuries Cost? A Systematic Review of the Economic Evaluations of Needlestick and Sharps Injuries Among Healthcare Personnel. Infect Control Hosp Epidemiol. 2016;37(6):635-646. doi:10.1017/ice.2016.48