

VOCO Perfect Protect®

DESENSITISING GEL FOR USE IN TOOTH WHITENING

Extension of the Perfect Bleach family

Many people would like whiter teeth without the discolourations, which can develop if they drink coffee or smoke, for example. VOCO has already been offering the tried-and tested, very high-quality products of the Perfect Bleach range for tooth whitening for years. Home bleaching methods in particular are very popular, as they can be optimally integrated into users' everyday lives. But what happens if teeth become hypersensitive as a result of bleaching at home?

The perfect complement to home bleaching!

VOCO Perfect Protect is a desensitising gel for use at home, which is tailored to the conditions of home bleaching. It contains 3% potassium nitrate (\triangleq 11,600 ppm potassium) and 0.11% fluoride (\triangleq 1,130 ppm fluoride), thus reliably contributing to the desensitisation of sensitive teeth. While the potassium nitrate aids with pain relief in hypersensitive teeth, the fluoride supports the sealing of open dentinal tubules and thus addresses the cause of the hypersensitivity.

The transparent and tasteless material is extremely simple and intuitive in its use: it is applied for 15–60 minutes using the same tray as for home bleaching and can be applied both before and after tooth whitening.

Versatile:

VOCO Perfect Protect is compatible with all home bleaching materials available on the market – meaning you do not need to be a Perfect Bleach customer to reap all the benefits of VOCO Perfect Protect!

The usual high efficiency:

VOCO Perfect Protect underscores the usual great performance of VOCO products once again. Measurements by Dr. Dimitrios Dionysopoulos from the department of Operative Dentistry at Aristotle University of Thessaloniki (Greece) showed that VOCO Perfect Protect performed best in the testing for the sealing of open dentinal tubules among the desensitising products intended for patient use. This means that VOCO Perfect Protect with around 80% sealed dentinal tubules (after use) is most effective against hypersensitivity caused by home bleaching (s. chart on the right).

Indication

Sealing of exposed dentine


Advantages

- Reliable desensitisation
- Short and variable wearing time
- Neutral taste
- Pleasant consistency
- Easy filling of the tray

Presentation

REF 1019 Syringe 2 x 2 ml

Quantity of sealed dentinal tubules

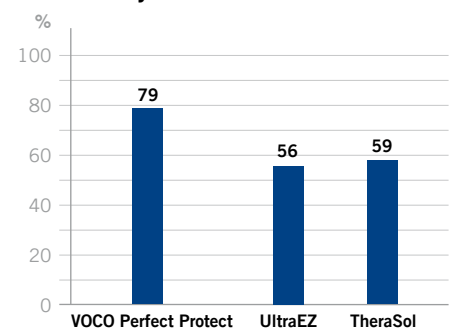

Source: Dr. Dimitrios Dionysopoulos, Department of Operative Dentistry, Aristotle University of Thessaloniki / Greece

Application



Clean teeth and fill material into the tray



Insert tray and remove excess material



Wear tray 15 – 60 minutes



Clean tray and teeth

VOCO GmbH
Anton-Flettner-Straße 1-3
27472 Cuxhaven
Germany

Tel.: +49 (0) 4721-719-0
Fax: +49 (0) 4721-719-140
info@voco.com
www.voco.dental

